

55" HD Medical Display




When Clarity Matters Most

At FSN, we understand that although space is limited in the OR, quality must not be compromised. Our complete medical monitor selection offers functionality, reliability, and all major inputs: DVI, VGA, Component, Composite, S-video, SOG, HDSDI and fiber.

- LED technology
- High definition picture
- Designed for medical use
- Rugged, IPX1 liquid resistant housing
- Withstands repeated cleaning
- Mounts securely to booms or custom yokes.
- Firmware can be upgraded and customized




Medical monitors from FSN are fully compatible with today's medical video providers.

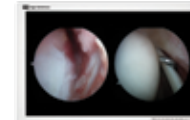


Multi-Display Modes

Picture-in-picture, picture-by-picture, zoom, swap images, pan, OSD lock and freeze functions can be controlled from the keypad located on the front of each display.

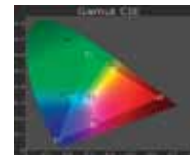

Accurate Color Reproduction

FSN monitors have the ability to clamp down the visible color spectrum, ensuring that doctors see exactly what they expect to see. Each display is calibrated to well known surgical color preferences. In addition, fine tune adjustments can be made at any time.

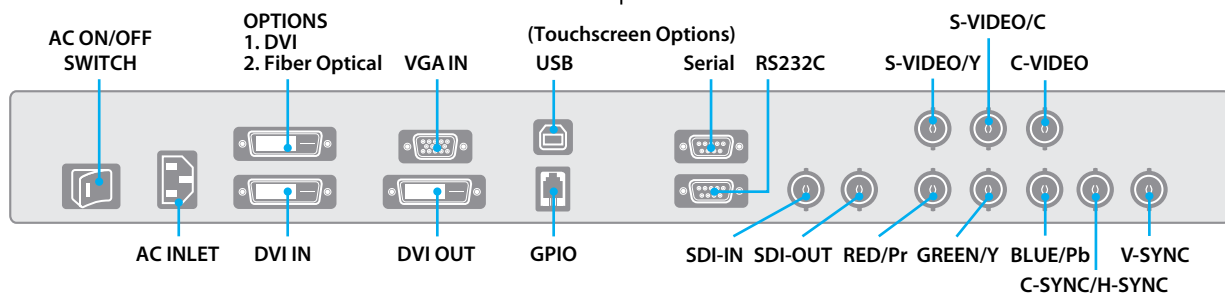

FSN displays have the connections you need. All major inputs are accepted via the back side connection panel. There are even connections for optional touchscreen capability.

Specifications





General Features

Item	Description
Panel	55" (LG Display) LED
Resolution	1920 x 1080 pixel
Input Signal	1 x DVI 1 1 x DVI 2 (optical fiber optional) 1 x D-SUB 1 x BNC (SDI) 1 x BNC (CVBS) 2 x BNC (SVHS Y/C) 5 x BNC (Component Y/G, Pb/B, Pr/R, H/CS, VS Input)
TOUCH (optional)	(Note: Depth and weight will increase)
GPIO port	Swap, PIP / PBP1 / PBP2 Select, Record indicator
Output Signal	1 x DVI & 1 x BNC (SDI)
Power Supply	AC 100~240V, 50/60 Hz, 3A Max
Control Key	Power, Menu, PIP, Up, Down, Plus, Minus, Input
Unit Dimension	1293(W) x 777(H) x 86.3(D) mm 50.9(W) x 30.59(H) x 3.4(D) inch
Average Brightness	500 (5 points avg.) cd/m2
Response Time (G-to-G)	6 ms (avg.)
Viewing Angle (CR>10)	R/L 178, U/D 178
Weight	37 kg, 81.57 lbs
Contrast Ratio	1400 : 1 typical
Pixel Pitch	0.630 x 0.630

Signal Timing

Input Signal	Type	Description
VGA	Type	Analog RGB
	Connector	DSUB-15
	Level	0.7Vp-p ±5%
	Polarity	Positive
	Impedance	75 ohm ±5%
	Horizontal frequency	30 ~ 93KHz
	Vertical frequency	50 ~ 85Hz
	Applicable maximum pixel frequency	170MHz
	Type	Separate H/V sync, Composite H/V sync, Sync On Green Automatic synchronization without external switch
	Terminal resistance	more than 2KΩ
DVI	Type	Digital RGB
	Signal Format	TMDS single link
	Connector	DVI-D
Optical DVI	Type	Digital RGB
	Signal Format	Fiber Optic Transmission
CVBS	Type	C-Video
	Color system	NTSC, PAL
	Signal Format	Composite Video Signal
SVHS Y/C	Type	S-Video
	Color system	NTSC, PAL
	Signal Format	Y/C Separated Video Signal
SDI		1080i (60 / 59.94 / 50)
	SMPTE-274M	1080p (30 / 29.97 / 25 / 24 / 24sF / 23.98 / 23.98sF)
	SMPTE-296M	720p (60 / 59.94 / 50)
	SMPTE-260M	1035i (60 / 59.94)
	SMPTE-125M	480i (59.94)
	ITU-R BT.656	576i (50)

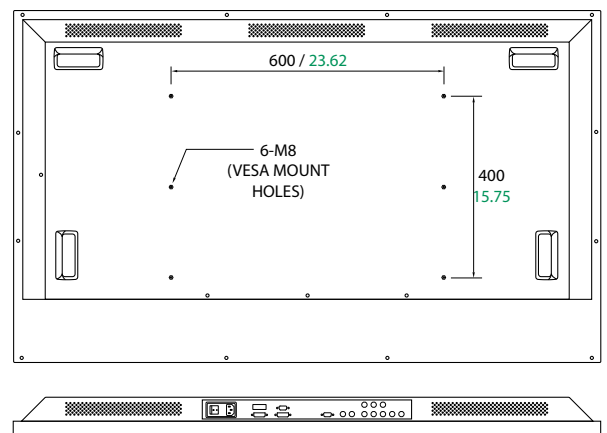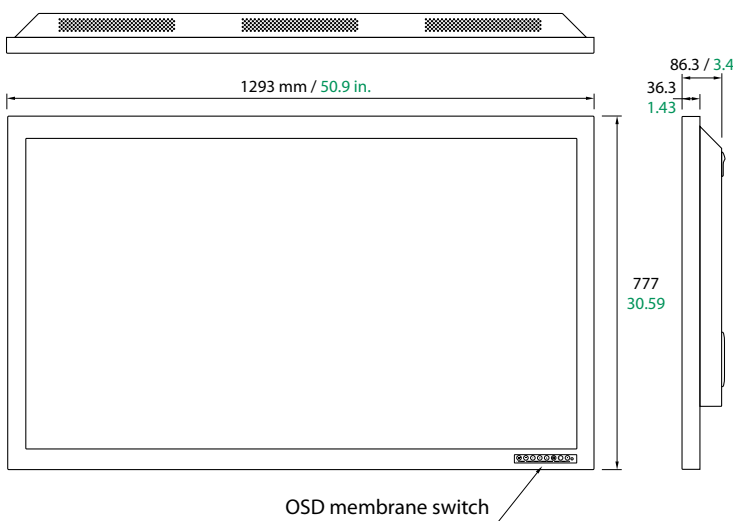

Specifications are subject to change with or without notice. Doc. # FSN1932 Rev. 10/12

www.fsnmed.com

Foreseon Custom Displays
2210 E. Winston Road
Anaheim, CA 92806
Tel: 714-300-0540
Fax: 714-300-0546

1 Bridge Plaza, Ste. 275
Fort Lee, NJ 07024
Tel: 201-849-4495
Fax: 201-490-1080

1800 Pembroke Dr., Ste. 300
Orlando, FL 32810
Tel: 407-667-3586
Fax: 407-667-4799

Benzstr.9
61352 Bad Homburg, Germany
Tel: +49(0)6172-185310-15
Fax: +49(0)6172-185310-11

Unit 2 Kingsmill Business Park
Chapel Mill Road
Kingston upon Thames, Surrey KT1 3GZ
Tel: 44 (0) 208 546 1047
Fax: 44 (0) 208 546 3931

Gang Nam Main Tower #801
1357-66 Seo-cho Dong
Seoul, 137-070 Korea
Tel: 82-2-521-5296
Fax: 82-2-521-5207

59-9 Jang-Dong, Yuseong-Gu
Daejeon City, Korea, 305-343
Tel: 82-42-360-8000
Fax: 82-42-360-8005